

Tables & Breadboards

Adjustable Breadboard Feet

These Breadboard feet can be used to stabilize any of our aluminum or stainless steel breadboards. They provide ±3mm (±0.120in) of vertical height adjustment to allow the breadboard to be leveled using a novel pivot design that compensates for uneven work surfaces. One foot is mounted under each corner to help reduce vibration. They are sold individually and are supplied with all necessary fixings.



Features


- ±3mm (±0.120") Vertical Height Adjustment for Easy Leveling
- Pivoting Foot for Uneven Surfaces
- Helps to Reduce Shock and Vibration

Breadboard Mounting Foot

ITEM#	METRIC ITEM#	\$	£	€	RMB	DESCRIPTION
BMF1	BMF1/M	\$ 6.12	£ 3.90	€ 5,70	¥ 58.40	Breadboard Mounting Foot, 1 Piece Only

Breadboard Lifting Handles

These handles enable breadboards to be lifted and transported easily. The glass fiber-reinforced polymer handle is safety rated to carry 250kg (550lbs). Each handle has two mounting holes that are used to bolt the handles to the breadboard and is supplied with all fixing screws. For larger breadboards, one handle at each corner is recommended.

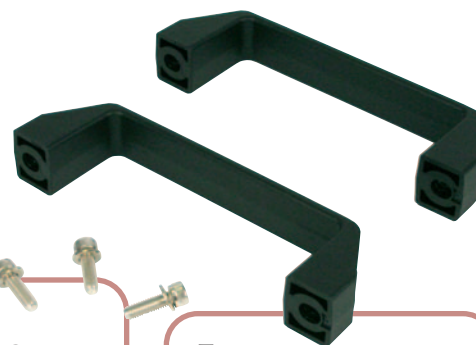

Specifications

- Each Handle Safety Rated to Carry 250kg (550lbs)
- Material: Glass Fiber-Reinforced Polyamide (PA) Technopolymer
- Resistant to Solvents, Oils, Greases, and Other Chemical Agents
- Color: RAL 9005 Black, Matte Finish
- Distance Between Mounting Holes: 150mm (6")

Features

- Easily Lift Breadboards From ScienceDesks and Work Frames
- Supplied With All Necessary Fixings
- Compatible With All Our Breadboards, Both Imperial and Metric




Breadboard Lifting Handles

ITEM#	METRIC ITEM#	\$	£	€	RMB	DESCRIPTION
BBH1	BBH1/M	\$ 16.95	£ 10.70	€ 15,80	¥ 161.90	Breadboard Lifting Handles (Set of 2)

Sorbothane Sheets

This pliable Sorbothane Sheet is ideal for isolating mechanical noise or vibrations. This substrate absorbs up to 95.5% of impact energy and then reforms to its original shape. The 1/4" or 1/8" thick sheets can be placed under your breadboard in single or stacked layers to counter the effects of uneven work surfaces or to isolate your breadboard from vibrations. For objects over 5lbs, the 1/4" thick sheets are recommended. The Sorbothane sheet may be hand-cut, die-cut, or kiss-cut to suit your applications.

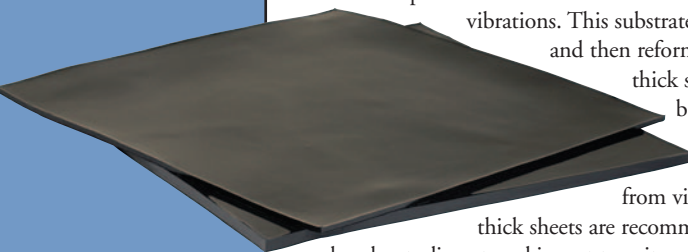

Features

- 12" x 12" Sheet Size With 1/4" or 1/8" Thickness
- Easy to Cut and Shape for Multiple Uses
- Hardness 70 Durometer (Type A)
- Absorbs Up to 95.5% of Impact Energy
- Can be Used as Vibration Isolation Material
- Ideal for Use With Our MB Series Aluminum Breadboards

Sorbothane Sheets

ITEM	\$	£	€	RMB	DESCRIPTION
SB12A	\$ 25.00	£ 15.80	€ 23,30	¥ 238.80	Sorbothane Sheet, 12" x 12" x 1/8", 70 Durometer
SB12B	\$ 49.00	£ 30.90	€ 45,60	¥ 468.00	Sorbothane Sheet, 12" x 12" x 1/4", 70 Durometer

Visit www.thorlabs.com For Mechanical Drawings and Our New Solid Models

